

Microfil® Filtration System

An easy to use, fast and accurate system for the routine microbiological analysis of your beverages and drinking water

The Microfil® system is a simple, low-cost and reliable solution for routine testing of raw materials, beverages and water for microbiological contamination. These tests are an important part of your quality assurance program. The standard testing method used in laboratories around the world for microbial analysis of fluids is membrane filtration. However, this method is technically demanding and requires time-consuming preparation, assembly and sterilization of the test apparatus.

With Merck Millipore's Microfil® filtration system, you will eliminate time-consuming steps, equipment set-up, waste, and autoclaving while improving the productivity of your lab. This cost-effective system uses presterilized, ready-to-use funnels and membranes on a filtration support.



Benefits

- No Preparation
- No Washing
- No Autoclaving
- No Clamps
- No Breakage
- No Filter Wrinkling
- No By-pass

Minimize Clean-Up, Assembly and Autoclaving

Ready-to-use, ultra-thin 100 or 250 mL Microfil® funnels eliminate the need for washing and sterilizing after each test and make equipment handling easy. Packaged in convenient stacks of 25 or 30, the funnels are removed as needed from the custom designed funnel dispenser. Such vertical stacking of the funnels also optimizes lab space.

A range of pre-sterilized 47 mm diameter gridded membranes are available with the Microfil® funnels. Each membrane is individually sealed in pleated band (EZ-Pak® membranes) or in peel back envelopes (S-Pak™ membranes). Polycarbonate membranes are also available.

Easy-to-Use with No Clamps and No O-Rings

The push-fit design of the Microfil® funnels seals tightly to manifold supports without the need for clamps or O-rings. This ensures leak-free operation and uniform microorganism recoveries.

Membrane Lifting Device

The membrane lifting device automatically breaks residual vacuum and lifts the membrane for easy removal with forceps. This aids in removing the membrane from the support after filtration, also reducing user-related contamination during the membrane transfer process.

Consistent Recoveries

Unlike with conventional filter holders, incomplete sealing at the membrane/funnel interface is eliminated resulting in uniform recovery of organisms. The smooth hydrophobic funnel surface repels sample residues and microorganisms to ensure that any microorganisms in your sample are collected on the membrane and not lost on the funnels walls. Pleating and distortion of the wet membrane due to expansion is avoided by the design of the support. As a result, uniform contact between the membrane and the medium is achieved.



Faster Filtration Rates

Due to the membrane support, filtration rates with the Microfil® system are faster when compared to traditional filtration equipment.

Regulatory Requirements

The Microfil® method has been established with reference to International Standards allowing biological analysis to be conducted under optimum conditions and conforms to International Standards for drinking water and mineral water.

ISO 7704 – Water Quality – Evaluation of membrane filters used for microbiological analysis – 1st edition 1985

Standard Methods for the Examination of Water and Waste Water 21st edition 2005

World Health Organization – Guidelines for drinking water quality, Geneva 1993


ISO 11731, Water quality: Detection and enumeration of *Legionella* spp, 1998 and 2004 (part 2)

Microfil® Method Advantages

- Use of 47 mm gridded presterilized membranes
- 100 mL and 250 mL funnels for drinking water and mineral water analysis
- Hydrophobic funnel surfaces for improved recovery
- Smooth interior funnel surfaces with no scratches
- Lip seal to prevent membrane by-pass
- Food-contact-polymer funnels, free from agents that may inhibit bacterial growth
- Disposable funnels eliminate decontamination by flaming or use of alcohol
- Recyclable plastic funnel material that allows for greater environmental protection
- 0.45 µm membranes are certified in accordance with Standard Methods for coliform analysis

Culture Media

You will find the full range of Merck Millipore culture media (including broth and dehydrated media) on our website at www.merckmillipore.com.

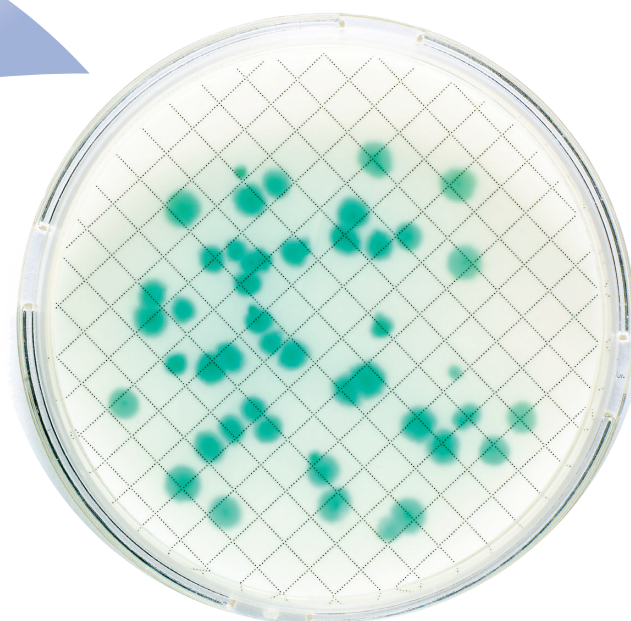

Processing an aqueous sample using the Microfil® system can be done in 7 simple steps:

Step 1

Sanitize the Microfil® support prior to processing each sample. Recommended methods include the use of alcohol or a quick flaming of the Microfil® support steel surface.

Step 2

Place either an EZ-Pak®, S-Pak™ or a polycarbonate membrane filter on the Microfil® support.



Step 3

Dispense a Microfil® funnel using the specially designed dispenser.



Step 4

Place the Microfil® funnel on the support and push down to fix it firmly in place.



Step 5

Pour the sample into the funnel and filter by applying vacuum.



The sealing lip of the Microfil® funnel avoids the risk of membrane by-pass by microorganisms. The clipping (1) holds the funnel firmly on the filtration support. After filtration, the funnel is removed and the membrane is lifted together with the membrane support (2) by simply pressing the lever (3). The membrane can now be conveniently removed from the support using forceps.

Step 6

Remove the funnel. Simply press the lever to lift the membrane off the support for easy removal with forceps. Residual vacuum is broken at the same time.



Step 7

Place the membrane filter into a Petri dish containing solid or liquid medium and incubate.



Ordering Information

150 Funnels and 150 S-Pak™ Membranes

Each box contains 6 x 25 (100 mL) or 5 x 30 (250 mL) sterilized Microfil® funnels and 150 individually packed, sterilized 47 mm diameter gridded S-Pak™ membranes, available with the following pore sizes and colors:

| Membrane Description | 100 mL Funnels | 250 mL Funnels |
|------------------------|----------------|----------------|
| 0.2 µm white, gridded | MIGSWG100 | MIGSWG250 |
| 0.45 µm white, gridded | MIHAWG100 | MIHAWG250 |
| 0.45 µm white, plain | MIHVWP100 | N/A |
| 0.45 µm black, gridded | MIHABG100 | MIHABG250 |
| 0.7 µm white, gridded | MIHCWG100 | MIHCWG250 |
| 0.8 µm white, gridded | MIAAWG100 | MIAAWG250 |
| 1.2 µm white, gridded | MIRAWG100 | MIRAWG250 |

150 Funnels and 150 Isopore Polycarbonate Membranes

| Membrane Description | 100 mL Funnels | 250 mL Funnels |
|--------------------------|----------------|----------------|
| 0.4 µm; colorless, plain | N/A | MIHTTP250 |

Additional membranes are available - 0.4 µm; colorless, plain - 150 units - item MIHTTPMNE

150 Funnels and 150 EZ-Pak® Membranes

Each box contains 6 x 25 (100 mL) or 5 x 30 (250 mL) sterilized Microfil® funnels and a band of 150 individually packed, sterilized 47 mm diameter gridded EZ-Pak® membranes, available with the following pore sizes and colors:

| Membrane Description | 100 mL Funnels | 250 mL Funnels |
|------------------------|----------------|----------------|
| 0.2 µm white gridded | MZGSWG101 | N/A |
| 0.45 µm white, gridded | MZHAWG101 | MZHAWG251 |
| 0.45 µm black, gridded | MZHABG101 | MZHABG251 |
| 0.8 µm white, gridded | MZAAWG101 | MZAAWG251 |
| 0.8 µm black, gridded | MZAABG101 | N/A |

Equipment

| Description (1/Pk) | Catalogue No. |
|--|---------------|
| 3-Position Vacuum Manifold, complete with Microfil® supports | MIAC03P01 |
| 1-Position Vacuum Manifold, complete with Microfil® supports | MIAC01P01 |
| Microfil® Support with integrated valve 1 1/4" NPTMI | MIAC01401 |
| Microfil® Funnel Dispenser, for use with 100 mL funnels | MIACFD101 |
| Microfil® Funnel Dispenser, for use with 250 mL funnels | MIACFD201 |
| EZ-Pak® dispenser | EZDISP001 |
| EZ-Stream™ Vacuum Pump for Liquid Transfer | EZSTREAM1 |
| Milliflex® plus vacuum pump base | MXPPUMP01 |
| Microfil® pump head for Merck Millipore vacuum pump | MCLHEAD01 |

Petri Dishes

47 mm Petri Dishes, packed in sleeves of 25 dishes. Available with or without cellulose pad.

| Description (1/Pk) | Qty/Pk | Catalogue No. |
|--|------------------------|---------------|
| Petri Dish without pad | 150 | PD2004700 |
| Petri Dish without pad | 600 | PD2004705 |
| Petri Dish with pad | 150 | PD2004750 |
| Petri Dish with pad | 600 | PD2004755 |
| 10 canisters of 100 absorbent pads | 1000 | AP1004550 |
| 1 pad dispenser and 2 pad canisters | 200 pads & 1 dispenser | AP1004551 |
| 0.45 µm with MCE membrane absorbent pads in canister | 600 | HAWG6745P |

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

For more information on our products:
www.merckmillipore.com/biomonitoring

Find contact information for your country at:
www.merckmillipore.com/offices

For Technical Service, please visit:
www.merckmillipore.com/techservice



Merck KGaA
Frankfurter Straße 250
64293 Darmstadt, Germany
e-mail: mibio@merckgroup.com